 $\oplus$ 

97050847-1 04/2015

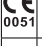
**GB** (translation from the original in Italian)

**Classification**  
Class Ila Medical device according to directive 93/42/EEC and subsequent additions and amendments. This device has been manufactured in compliance with IEC 60601-1 and ISO 14457.

**Intended Use**  
Medical device, or a medical device accessory for professional use by Dentists, assisted by authorised personnel. The device is intended for use on dental units complying with the 93/42/EEC Directive and subsequent additions and amendments complying with the IEC 60601-1 standard (III ed.) and installed in locations with an electrical plant system that complies with the IEC 60364-7-710.

**Symbols**

	Symbol to call the attention to other information found in the equipment User Manual.
	Symbol corresponding to "TYPE B PART APPLIED" according to IEC 60601-1. Indicates a protection grade against direct and indirect contacts
	Part that may be sterilised in an autoclave
	Mark of Compliance with the 93/42/EEC Directive "MEDICAL DEVICES" and subsequent additions and amendments.
	Symbol for disposal pursuant to the 2002/95/EC, 2002/96/EC, 2003/108/EC Directives.

**Precautions**  
•The device is not suitable for use with blends of inflammable anaesthetics and oxygen or nitrous oxide.  
•Use of the device must respect the indications contained the Instruction Manual available on the website [www.sternweber.com](http://www.sternweber.com).  
•The instrument is provided in a non-sterile state. Before use, see to its suitable cleaning, disinfection and/or sterilisation (as applicable).  
•Use of electrically-powered devices may interfere with the operation of active implantable devices such as pacemakers or other active devices. In case of doubt regarding the treatment of patients with such devices, consult a cardiology expert or another competent medical institute.  
•During use of the device dust and fragments from the patient's oral cavity or from the same device may be projected into the environment (organic and inorganic dust and metal dust, fragments of the device or its tips and potentially infected biological material).  
•Protect the patient, when possible, using a dental dam.  
•Instruct the patient to breathe through the nose in cases where the dam is not applicable.  
•Professional personnel must put on suitable personal protection equipment.  
•In case of visible damage or of the emission of unusual noises and/or vibrations or if overheating is noticed, do not use the device and call authorised technical service assistance.  
•Before each use, verify the proper locking of the burrs. If the burr does not fit correctly, do NOT use the tool.  
•Do not use worn or damaged burrs. Handle the burrs carefully wearing protective gloves.  
•The burr lock button overheats if pressed while the turbine is operating.  
•Be very careful not to inadvertently press it inside the patient's mouth.  
•Use only burrs (or other similar tools) with CE Mark of compliance with Directive 93/42/EEC and subsequent additions and amendments, made with biocompatible materials complying with ISO 10993-1 and whose dimensions comply with EN ISO 1797-1 standards.  
•Use only original spare parts and accessories, or those authorised by the manufacturer.  
•Make sure that the surgical area cooling flow is present during operation. In the absence of cooling, suspend use of the device and contact the authorized technical support service.  
•Do not breathe into beam or view directly with optical instruments.

**Assembling the burrs**  
Simply press the middle of the button on the head to enable the burrs to be inserted and removed. Releasing the button locks the burr. Check carefully that the burr is fully inserted and locked, pulling it with two fingers. Do not press the burr release mechanism if the turbine is not completely still in order to avoid damaging the burr locking mechanism.

**Periodic Maintenance**  
Use of the device must be carried out with respect for the indications found in the specific Instruction Manual.  
The turbine must be lubricated at least once a day before each time it is sterilized in the autoclave with the special "DAILY OIL" spray using the appropriate adapters. Follow the instructions on the can, making sure to apply the spray in two short jets. After lubrication, replace the turbine on the cord and activate it in order to expel the excess lubricant (do this with a burr inserted in order to avoid mechanical damage to the instrument), then dry the outside of the device with gauze or cotton. Regular use of the lubricant spray ensures the best operating conditions for the instrument. At least once a week, lubricate the "O-Rings" located on the quick coupling, using S1 silicone lubricant. Wearing disposable gloves, apply a layer of grease on the fingertips and lubricate using the fingers. If liquids that potentially leave incrustations are used to supply the spray, it is recommended that the circuit be flushed well with water at the end of each use. The manufacturer requires an annual check or overhaul of the device by authorised service assistance personnel.

**Cleaning, Disinfection And Sterilisation**  
Before disinfection, remove the burr from the turbine and the turbine from the cord. To maintain the level of hygiene, after each use and within a short time, clean, disinfect and sterilize the turbine. Clean the spray nozzles spray both by using the special steel wire supplied and by blowing in air using the syringe. For cleaning and disinfection of the outside, use gauze or cotton soaked in ethyl alcohol 70% v/v. Sterilise in a steam autoclave max 135° C 220 kPa (2.2 bar) for 5 minutes (250 cycles tested – EN ISO 14457). Before each cycle in the autoclave, it is essential to lubricate the turbine using the "DAILY OIL" lubricant spray, applying two short jets. Reassemble the turbine on the cord and activate it in order to expel the excess lubricant. Dry the outside of the turbine with gauze or cotton and remove the quick coupling, if present, prior to sterilization.

•Do NOT sterilize the quick coupling.  
•Do not use ultrasound cleaners.  
•Do not immerse the turbine in disinfectant or sterilizing liquids.  
•Never sterilise the handpiece in a dry heat steriliser.  
•Do not leave the turbine in the autoclave at the end of the cycle.  
•Periodically check the autoclave according to the requirements of the manufacturer.

Sterilization temperatures above the limit may damage the device.

**Repairs and overhauls**  
Installation, overhaul, calibration and repair of the device must be performed exclusively by technical service assistance personnel authorised by the manufacturer.

**Intended Use**  
The method of application of the device is intermittent. The operational cycle indicated respects the following times:

Type of Instrument	Operation (minutes)	Pause (minutes)
Silent Power EVO 2	5	5
Silent Power EVO 2L	5	5
Silent Power EVO 4L / 4L-K	5	5
Silent Power EVO Miniature	5	5
Silent Power EVO Miniature K	5	5

**Environmental conditions for use:**  
- Ambient Temperature 10 - 40° C  
- Relative Humidity 30 - 75 %  
- Atmospheric Pressure 700 - 1060 hPa (700 - 1060 mBar)

**Environmental conditions permitted for transport and warehousing**  
- Ambient Temperature between -20 - +70° C  
- Relative Humidity 10 - 100 %  
- Atmospheric Pressure 500 - 1060 hPa (500 - 1060 mBar)

**Scraping**  
Scraping the device must be carried out with respect for the laws in force and according to the individual national laws. The materials used for the manufacture are not harmful for human beings or animals that come into contact with or that are exposed to them.

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